

Delta® Ultra-Light 1000, Full-Electric Bed

HCPCS:



FIGURE A



FIGURE B

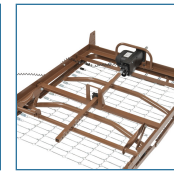


FIGURE C



FIGURE D



FIGURE E

- The Delta 1000 is truly universal because the bed frame can be used with Drive's and most other manufacturers' bed ends, old and new
- The transition box mounted on the foot section is key to interchangeability (Figure A)
- The transition box allows you to change the rotation of the high/low shaft
- Headboard is taller than footboard to better fit home decor (Figure B)
- The wood grain panels are break and scratch resistant
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations (Figure C)
- Labels and color-coded springs assure accurate installation of side rails for patient safety (Figure D)
- Unique motor is completely self-contained to reduce weight and noise
- The motor assembly can be installed and removed with patient in bed
- Bed, side rails and mattress safety-tested by SGS, a nationally recognized testing facility Satisfies all applicable FDA entrapment guidelines
- The motor is UL approved
- In the event of a power failure, a 9 volt battery, located in the hand control, can lower the head and foot sections nine times No crank needed
- New and improved hand pendant has large, easy-to-use controls (Figure E)
- Bed ships in two cartons
- Channel frame construction provides superior strength and reduced weight
- Zinc-coated spring deck
- Extension kit available (model #15030EXTKIT) extends length from 80" to 84"
- Patent Pending
- Weight Capacity 450 lbs

PRODUCTS	Description
15033	Full Electric Ultra Light Plus,1/cs
15033BV-FR	Ultra Lt.Full Elec.Bedw/Full R
15033BV-HR	Ultra Lt.Full Elec.Bed,w/Half
15033BV-PKG	Ultra Lt.Full Elec.Bed Pkg,1pk
15033BV-PKG-1	Ultra Lt.Full Bed Pkg w/HalfRa
15033BV-PKG-11	Delta Full,FullRail,QuickEasyFoam,1ea
15033BV-PKG-12	Delta Full,HalfRail,QuickEasyFoam,1ea
15033BV-PKG-1-T	Ultra Lt.Full Bed Pkg w/BrownV
15033BV-PKG-2	Ultra Lt.Full Elec.BedPkg w/Fo
15033BV-PKG-4	Full Elect. Ultra Lt. Plus Bed,1/cs
15033BV-PKG-T	Ultra Lt.Full Elec.Bed Pkg,1pk

SPECIFICATIONS	Description
Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components
Weight Capacity	450 lbs Maximum weight limit includes weight of occupant
Carton Shipping Weight	182 lbs.
Weight	169 lbs.
ACCESSORIES	Description
15030HBL	Hospital Bed Linen Kit,1 EA,10/CS, CS
15030HBC	Bedding in a Box,1 EA,5/CS, CS
PARTS	Description
15044CSET	Set of 4 Casters,Homecare Bed, 1 set/ea, EA
15544H-F	Head&Footboard15030/15033Delta,1ea, EA
15033MHLN-S	Mounting Screws,15033MHLN,1 set/bag, BAG
15005P1018	Pin,SpringAsbly,Semi,Full,1/ea, EA
15033HC-G	HandControlFullUltraLight,1ea, EA
15033MHLN-MGBA	Hi Lo Motor with Gear Box,1ea, EA
15549	Foot Spring,15033,15235,Delta,1/ea, EA
15037MN	Motor,Head/Footspring,15033/15235,1/ea, EA
15033MBN-DELTA	Hi-LowMotorMountingBracketDelta,1ea, EA
15544HGB	Head Board Gear Box,15544,1ea, EA
15044P	Replacement BedEnd Panel,15044,1/ea, EA
15544FPF	Foot BoardPanelFrame,DarkBrown,15544,1ea, EA
15541	Head Spring,15030/15033,Delta Bed,1ea, EA
15544H	HeadBoard,15030/15033,Delta,Bed,1ea, EA
15544HPF	Head BoardPanelFrame,DarkBrown,15544,1ea, EA
15030-FRLINK	FrameLinkw/DeckBolts,15030,15005, 1/ea, EA
15033MHLNC-N	HiLoMotorConnector15033MHLN-MGBA,1ea, EA
15544GBMB	MountingBracketForTransBox15030,1ea, EA
15544FGB	Foot Board Gear Box,15544,1ea, EA
15544F	Foot Board,15030/15033,Delta Bed,1ea, EA