

Deluxe, Trigger Release Folding Walker with 5" Wheels

HCPCS: E0143,E0143



FIGURE A

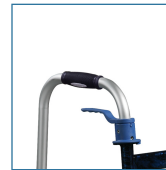


FIGURE B

- The Deluxe Trigger Release Folding Walker has a trigger release feature that allows the user to close the walker without removing their hands from the handle, making it ideal for individuals with limited finger dexterity. The lightweight aluminum frame ensures maximum strength and durability. The 5" wheels allow for navigation on many types of terrain, while the contoured hand grip and adjustable height make the walker comfortable to use.
- Durable, Composite Trigger Release resists cracking and breaking
- Sturdy aluminum construction ensures maximum strength while remaining lightweight
- Trigger release feature for opening and closing allows users to close the walker without removing their hands from the handle and is ideal for individuals with limited finger dexterity
- Durable, composite trigger release can be reached with hands of all sizes and resists cracking and breaking, ensuring safety and durability
- Ribbed, contoured handgrips relieve pressure on users hands
- Using the marked peg holes, easily adjust the height in 1" increments for a customized fit
- 5" wheels are ideal for maneuvering both indoors and outdoors
- Rear glide cap allows for walker to slide easily and smoothly over most surfaces
- Folds for easy transport and the sides are detachable for storing in tight spaces
- Quick and tool-free assembly

PRODUCTS	Description
10226-1	Walker,TriggerRls,Adlt,w/5"Whl,1/ea
10226-4	Walker,2 Trigger Rel.Adlt.w/5",4/cs
10227-4	Walker,2TriggerRel.Junr.w/5"Whls,4/cs
SPECIFICATIONS	Description
Warranty	Limited Lifetime
Material	10226-4: Aluminum
Handle Height	10226-1 and 10226-4: 32.5" - 39.5" 10227-4: 25.5" - 33"
Weight Capacity	350 lb
Carton Shipping Weight	10226-1 and 10226-4: 29 lb 10227-4: 26 lb
Color	10226-4: Silver
Weight	10226-1 and 10226-4: 6 lb 10227-4: 6 lb
Wheel Size	10226-1, 10226-4 and 10227-4: 5"
CA Prop 65	[WARNING: This product can expose you to chemicals that are known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov]
ACCESSORIES	Description
10200B	Walker Basket,1/ea, RTL, EA

ACCESSORIES	Description
RTL10131	Walker Caddy, 3/cv, CS
10121	Tennis Ball Glides, 1pr/bx, RTL, BX
10110	Glide Ski for Walker, 1pr./bx, RTL, BX
RTL10254BFL	Univ Mobility Bag, Med, Blue Flrl, 1EA, EA
RTL10254GR	Univ Mobility Bag, Med, Army Green, 1EA, EA
RTL10254TFL	Univ Mobility Bag, Med, Trop Flrl, 1EA, EA
RTL10255NV	Univ Mobility Bag, Small, Navy, 1EA, EA
RTL10255PK	Univ Mobility Bag, Small, Pink, 1EA, EA
RTL10256NV	Univ MobilityBag, Med Flat, Navy, 1EA, EA
RTL10256PK	Univ MobilityBag, Med Flat, Pink, 1EA, EA