

## Jazzy® Carbon<sup>1</sup>

HCPCS Code	K0899 Power Mobility Device, not coded by DME PDAC or does not meet criteria
Weight Capacity	300 lbs.
Turning Radius <sup>3</sup>	24"
Overall Width <sup>3</sup>	24"
Overall Length <sup>3</sup>	37"
Maximum Speed <sup>2</sup>	Up to 3.7 mph
Ground Clearance <sup>3</sup>	2"
Front Casters <sup>3</sup>	7"
Drive Wheels <sup>3</sup>	8.5" solid
Rear Casters <sup>3</sup>	N/A
Drivetrain	24v/250W brushless motors
Braking System	Regenerative and electromechanical
Suspension Type	Front
Standard Seat <sup>3</sup>	17" x 17" black foam
Seat-to-Floor Range <sup>3</sup>	20.25" with standard cushion 18.75" without cushion
Per Charge Range <sup>3,5</sup>	Up to 9.3 miles
Standard Electronics	35A Brushless Micon
Component Weights <sup>6</sup>	Base: 39.6 lbs. Battery: 4 lbs.
Total Weight <sup>4</sup>	43.6 lbs. with batteries
Battery Requirements	24V 12AH (288WH) lithium-ion
Battery Charger	2 amp, lithium-ion off-board
Warranty	Lifetime limited on frame; 1-year limited on electrical control/ lever system, motor/drive system, bearing sleeve and battery



FDA Class II  
Medical Device<sup>1</sup>

**PRIDE**  
MOBILITY

800-800-8586 | pridemobility.com



HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at time of publication; we reserve the right to alter specifications without prior notice.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

- <sup>1</sup> Jazzy® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.
- <sup>2</sup> Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).
- <sup>3</sup> Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- <sup>4</sup> Total weight includes base, seat and batteries.
- <sup>5</sup> Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- <sup>6</sup> Battery weight may vary based on manufacturer.

**Jazzy**  
**carbon**

**FOLDS FLAT  
FOR TRANSPORT**

**OPENS  
WITH  
ONE-HAND  
SNAP  
ACTION!**

**LIGHT  
ENOUGH TO  
PICK UP!**

**AMERICA'S #1 TRAVEL POWER CHAIR**

pridemobility.com |    

FDA Class II Medical Device<sup>1</sup>



# Jazzy<sup>®</sup> carbon

39.6  
LBS  
WEIGHT WITHOUT BATTERY

Padded seat back  
with storage pocket

Quick disconnect joystick  
with integrated USB charger

Mesh Cup holder

Ultra-lightweight,  
carbon fiber frame


Compact base easily  
folds and unfolds for  
easy transport in a  
few easy steps

Airline compliant  
lithium battery\*

Front suspension  
for a smooth ride

Under-seat  
storage

FDA Class II  
Medical Device<sup>1</sup>

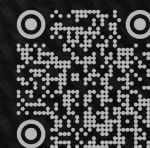

## UNFOLD, SNAP, AND GO!

This is our most portable Jazzy<sup>®</sup> yet. In just seconds, the Jazzy Carbon folds flat for easy transport, and unfolds again to get you moving. Ultra-light but exceptionally strong, this power wheelchair is your ultimate travel companion.

Weighing less than 40 lbs. (without battery), the Jazzy Carbon still provides top speeds up to 3.7 mph with a weight capacity of 300 pounds.

### FEATURES

- Quick disconnect joystick with built in USB charger
- 24V 12Ah (288WH) airline compliant lithium battery\*
- Front suspension for a smooth ride
- Ultra-lightweight, carbon fiber frame specifically designed for traveling
- Compact base easily folds and unfolds for transport in a few simple steps
- Padded seat back with storage pocket and soft, removable seat cushion for all-day comfort
- Free-wheel levers easily disengage the brake, allowing this lightweight power chair to be pushed like a manual chair
- Under-seat storage and adjustable lap belt
- Standard mesh cup holder

### AVAILABLE COLORS



Black

White

### ACCESSORIES

- Additional battery
- Attendant control bracket
- **NEW!** Battery bag
- Headrest
- **NEW!** Jazzy<sup>®</sup> Travel Kit (includes a bag to house the folded unit and a separate battery bag)
- Military patches
- Stealth<sup>®</sup> Simplicity cushion

**NEW!**



Battery Bag



Height: 28.3" | Width: 24" | Depth: 12.2"



The Jazzy<sup>®</sup> Travel Kit includes a bag to protect your folded Jazzy<sup>®</sup> Carbon and a separate bag for the battery.

Folds quickly and easily for transport or storage.

Front suspension provides a smooth ride.

<sup>\*</sup> Transport of this item on commercial aircraft is subject to federal regulations and airline policy. Please consult your airline prior to your flight to ensure your Jazzy<sup>®</sup> Carbon and its lithium ion battery are permitted.